


Clinical assessment and self-assessment bench Various, stable and calibrated

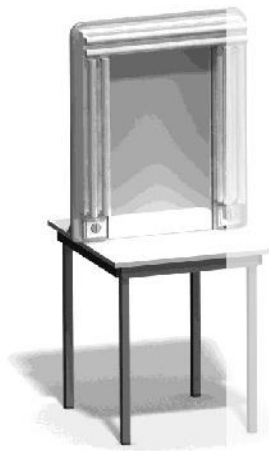

→ INTRODUCTION

The Evalux Bench helps to perform clinical assessment of face, low neck, hair and hands in repeatable conditions thanks to stable and calibrated lightings. Thanks to a mirror, it allows self-assessment of the panelist or assessment by a clinical expert according to same conditions such as positioning or illumination. This flexible system can be adapted to various user requirements. The various types of illuminants allow to reproduce in a reliable way some environmental conditions such as «Daylight», solar or artificial lightings for various incidences on the face (Lighting «through the top», symmetrical right/ left or the combination of both).

→ CHARACTERISTICS

LIGHTINGS :

Homogeneous, without shadows, on the whole face

- Perfect and complete stability of the illumination
- Long lifetime (50000h)
- Powerful lighting minimizing the influence of ambient lighting
- 3 tubes with 180 white super luminescent LEDs
- Possible illuminants : 6000K ° (daylight), 3100K ° (solar light) and others
- IRC 96 (Index of optimal depiction of colors)
- Continuous spectrum
- Possibility of 3 supplementary tubes for the addition of other

EVALUATION :

- Mirror 30X40 cm for self-assessment
- A mirror foldaway for the expert clinical assessments (Sliding door)
- Possibility of lighting from top (1 tube); symmetrical right/left (2 tubes) or combination of both lighting incidences (3 tubes)
- Space dedicated to a laptop

MAGNA ALIQUAM ERAT VOLUTPAT	
LOREM IPSUM DOLOR	Lorem ipsum Laoreet dolore magna
CONSECTETUER	Sit amet, consetetuer adipiscing Laoreet dolore magna aliquam erat volutpat
ADIPISCING ELIT	Sed diam nonummy Laoreet dolore magna aliquam Consetetuer adipiscing elit

Figure

